

# InstrumentLab 7.5

## Quick Start



[www.openwire.org](http://www.openwire.org)  
[www.mitov.com](http://www.mitov.com)

Copyright Boian Mitov 2004 - 2014

# Index


<a href="#">Installation.....</a>	<a href="#">3</a>
<a href="#">Where is InstrumentLab?.....</a>	<a href="#">3</a>
<a href="#">Creating application with angular gauge.....</a>	<a href="#">3</a>
<a href="#">Adding component elements to a gauge.....</a>	<a href="#">5</a>
<a href="#">Using the TSLCRealBuffer in C++ Builder and Visual C++.....</a>	<a href="#">8</a>
<a href="#">Deploying your 32 bit application with the IPP DLLs.....</a>	<a href="#">10</a>
<a href="#">Deploying your 64 bit application.....</a>	<a href="#">11</a>

# Installation

InstrumentLab comes with an installation program. Just start the installation by double-clicking on the Setup.exe file and follow the installation instructions.

## Where is InstrumentLab?


After the installation, start your Delphi or C++ Builder. Scroll the Component Palette, until you see the last three tabs:

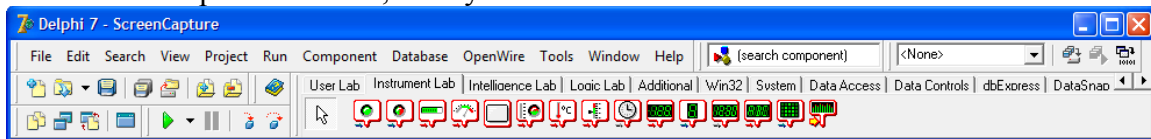

If the installation was successful, they should be named “Signal Lab”, “UserLab”, and “InstrumentLab”. On the SignalLab palette you will have only a subset of the SignalLab components. SignalLab is a separated product, and will not be shipped as full with InstrumentLab.


Only the following components of SignalLab will be available:

The following UserLab components will also be available:

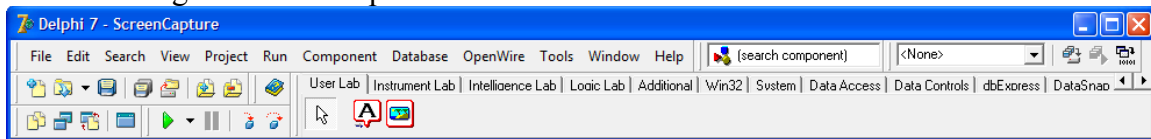


## Creating application with angular gauge

From the Delphi/C++Builder menu select | File | New | Application |.

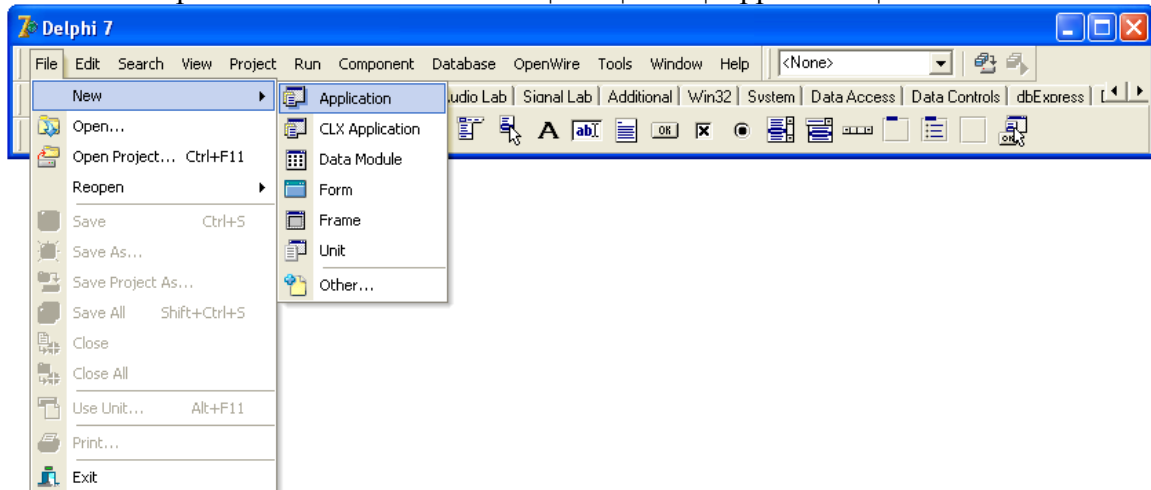

An empty form will appear on the screen.

From the “Component Palette” select the “Instrument Lab” tab:

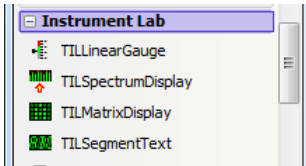


select and drop on the form the following component:



- TILAngularGauge

From the “Component Palette” select the “System” tab:

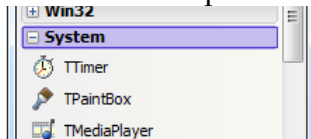


From the tab select and drop on the form a TTimer component.




- TTimer

Select the Timer1 component on the form:

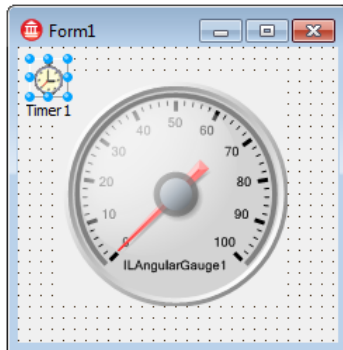

In the “Object Inspector” set the Interval to 100:

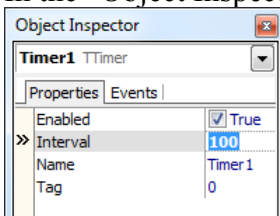

Double-click on the Timer1.

**If you are using Delphi add the highlighted code to the Timer1Timer event in the Unit1.pas file:**


```
procedure TForm1.Timer1Timer(Sender: TObject);  
begin  
    IAngularGauge1.Value := Random( 100 );  
end;
```

If you are using C++ Builder add the highlighted code to the Unit1.cpp file:

```
void __fastcall TForm1::Timer1Timer(TObject *Sender)
{
    IAngularGauge1->Value = rand() * 100 / RAND_MAX;
}
```

Compile and run the application.

You should see gauge showing the random values:

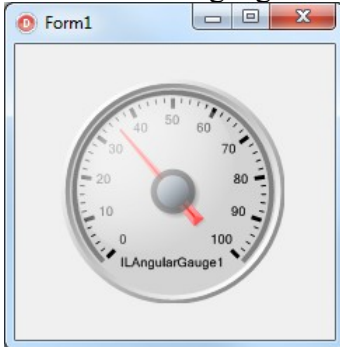


You have just learned how to use InstrumentLab gauges.


## Adding component elements to a gauge

Open the application already created in the “Creating application with angular gauge” chapter.


Double-click the Gauge1:

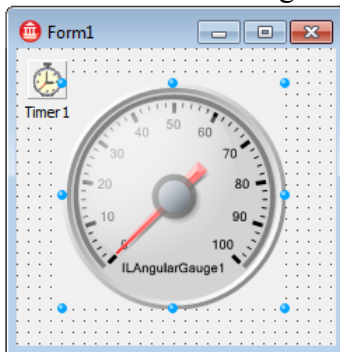

The “Component Editor” will appear:

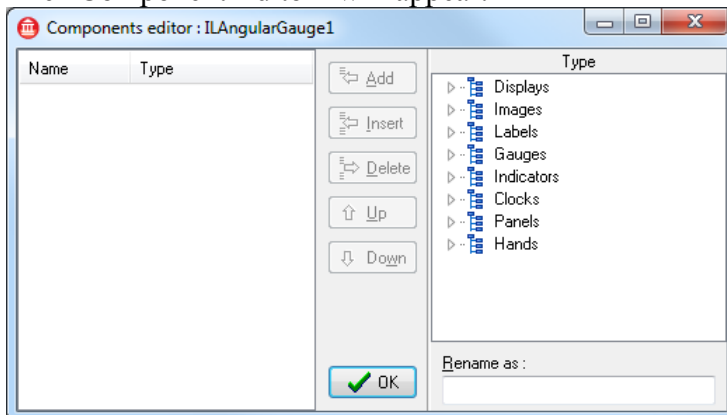

In the “Component Editor” add TILExternalGaugeHandElement:

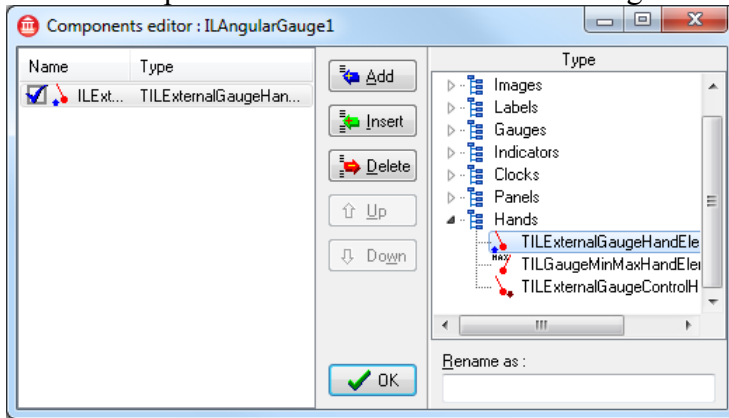

Click “OK”.

**If you are using Delphi add the highlighted code to the Timer1Timer event in the Unit1.pas file:**


```
procedure TForm1.Timer1Timer(Sender: TObject);
begin
  ILAngularGauge1.Value := Random( 100 );
  ILExternalGaugeHandElement1.Value := Random( 100 );
end;
```

**If you are using C++ Builder add the highlighted code to the Unit1.cpp file:**

```
void __fastcall TForm1::Timer1Timer(TObject *Sender)
{
  ILAngularGauge1->Value = rand() * 100 / RAND_MAX;
  ILExternalGaugeHandElement1->Value = rand() * 100 / RAND_MAX;
}
```


Compile and run the application.

You should see gauge showing the random values with both the hands:

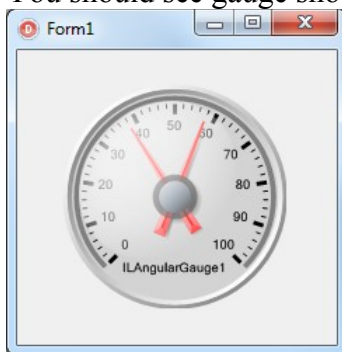


Stop the application.


Double-click the Gauge1:

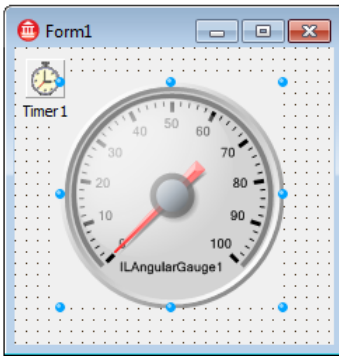


In the “Component Editor” add TILThermometerElement:

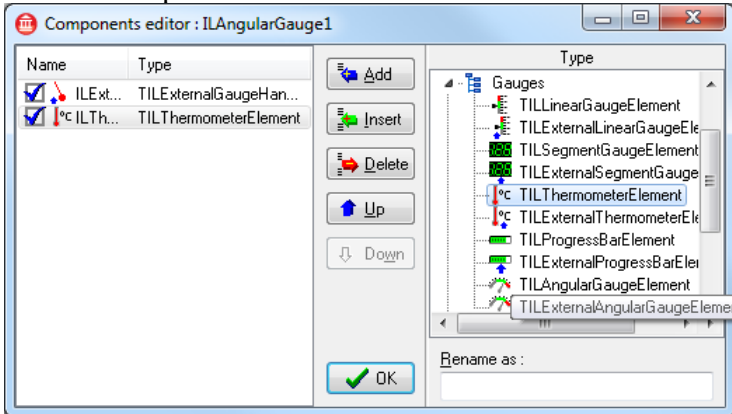

In the “Component Editor” add TILLevelDetectElement:

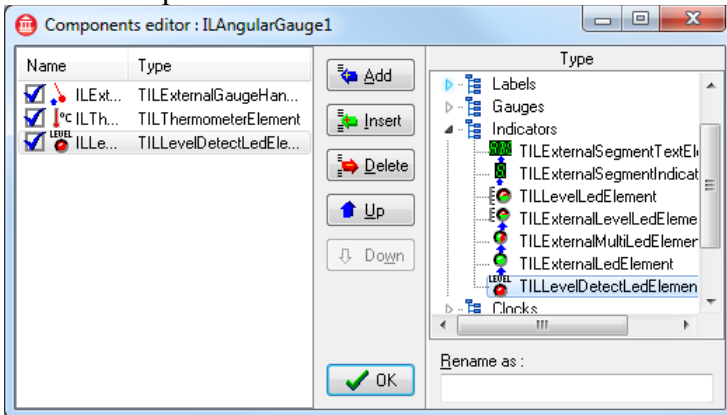


Click on the ILLevelDetectElement1 to show its properties in the “Object Inspector”.  
In the “Object Inspector”, set the Level property to 50:

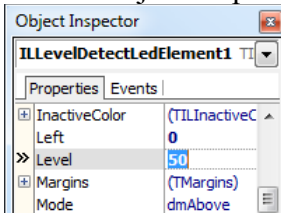

Compile and run the application.

You should see gauge showing the random values with both the hands, the LED will be On when the first hand is above 50 and the Thermometer will show the value of the first hand:

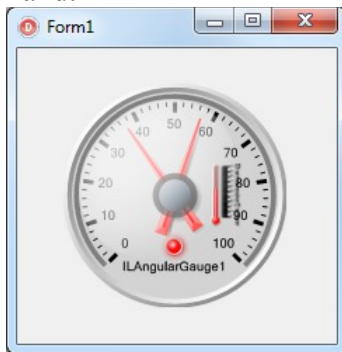

You have just learned how to create composite components in InstrumentLab.

## Using the TSLCRealBuffer in C++ Builder and Visual C++

The C++ Builder version of the library comes with a powerful data buffer class, called TSLCRealBuffer.

The TSLCRealBuffer is capable of performing basic math operations over the data as well as some basic signal processing functions. The data buffer also uses copy on write algorithm improving dramatically the application performance.

The TSLCRealBuffer is an essential part of the SignalLab generators and filters, but it can be used independently in your code.

You have seen already some examples of using TSLCRealBuffer in the previous chapters. Here we will go into a little bit more details about how TSLCRealBuffer can be used.

In order to use TSLCRealBuffer you must include SLCRealBuffer.h directly or indirectly (through another include file):

```
#include <SLCRealBuffer.h>
```

Once the file is included you can declare a buffer:

Here is how you can declare a 1024 samples buffer:

```
TSLCRealBuffer Buffer( 1024 );
```

Version 4.0 and up does not require the usage of data access objects. The data objects are now obsolete and have been removed from the library.

You can obtain the current size of a buffer by calling the GetSize method:

```
Int ASize = Buffer.GetSize(); // Obtains the size of the buffers
```

You can resize (change the size of) a buffer:

```
Buffer.Resize( 2048 ); // Changes the size to 2048
```

You can set all of the elements (samples) of the buffer to a value:

```
Buffer.Set( 30 ); // Sets all of the elements to 30.
```



You can access individual elements (samples) in the buffer:

```
Buffer [ 5 ] = 3.7; // Sets the fifth element to 3.7  
  
Double AValue = Buffer [ 5 ]; // Assigns the fifth element to a variable
```

You can obtain read, write or modify pointer to the buffer data:

```
const double *data = Buffer.Read() // Starts reading only  
double *data = Buffer.Write() // Starts writing only  
double *data = Buffer.Modify() // Starts reading and writing
```

Sometimes you need a very fast way of accessing the buffer items. In this case, you can obtain a direct pointer to the internal data buffer. The buffer is based on copy on write technology for high performance. The mechanism is encapsulated inside the buffer, so when working with individual items you don't have to worry about it. If you want to access the internal buffer for speed however, you will have to specify up front if you are planning to modify the data or just to read it. The TSLCRealBuffer has 3 methods for accessing the data Read(), Write(), and Modify (). Read() will return a constant pointer to the data. You should use this method when you don't intend to modify the data and just need to read it. If you want to create new data from scratch and don't intend to preserve the existing buffer data, use Write(). If you need to modify the data you should use Modify (). Modify () returns a non constant pointer to the data, but often works slower than Read() or Write(). Here are some examples:

```
const double *pcData = Buffer.Read(); // read only data pointer  
  
double Value = *pcData; // OK!  
*pcData = 3.5; // Wrong!  
  
double *pData = Buffer.Write(); // generic data pointer  
  
double Value = *pData; // OK!  
*pData = 3.5; // OK!
```

You can assign one buffer to another:

```
Buffer1 = Buffer2;
```

You can do basic buffer arithmetic:

```
TSLCRealBuffer Buffer1( 1024 );  
TSLCRealBuffer Buffer2( 1024 );  
TSLCRealBuffer Buffer3( 1024 );  
  
Buffer1.Set( 20.5 );  
Buffer2.Set( 5 );  
  
Buffer3 = Buffer1 + Buffer2;  
Buffer3 = Buffer1 - Buffer2;  
Buffer3 = Buffer1 * Buffer2;  
Buffer3 = Buffer1 / Buffer2;
```

In this example the elements of the Buffer3 will be result of the operation ( +,-,\* or / ) between the corresponding elements of Buffer1 and Buffer2.

You can add, subtract, multiply or divide by constant:

```
// Adds 4.5 to each element of the buffer
Buffer1 = Buffer2 + 4.5;

// Subtracts 4.5 to each element of the buffer
Buffer1 = Buffer2 - 4.5;

// Multiplies the elements by 4.5
Buffer1 = Buffer2 * 4.5;

// Divides the elements by 4.5
Buffer1 = Buffer2 / 4.5;
```

You can do “in place” operations as well:

```
Buffer1 += Buffer2;
Buffer1 += 4.5;

Buffer1 -= Buffer2;
Buffer1 -= 4.5;

Buffer1 *= Buffer2;
Buffer1 *= 4.5;

Buffer1 /= Buffer2;
Buffer1 /= 4.5;
```

Those are just some of the basic buffer operations provided by SignalLab.

If you are planning to use some of the more advanced features of TSLCRealBuffer please refer to the online help.

SignalLab also provides TSLCComplexBuffer and TSLCIntegerBuffer. They work similar to the TSLCRealBuffer but are intended to be used with Complex and Integer data. For more information on TSLCComplexBuffer and TSLCIntegerBuffer please refer to the online help.

## Deploying your 32 bit application with the IPP DLLs

The compiled applications can be deployed to the target system by simply copying the executable. The application will work, however the performance can be improved by also copying the Intel IPP DLLs provided with the library.

The DLLs are under the [install path]\LabPacks\IppDLL\Win32 directory( [install path] is the location where the library was installed).

In 32 bit Windows to deploy IPP, copy the files to the [Windows]\System32 directory on the target system.

In 64 bit Windows to deploy IPP, copy the files to the [Windows]\SysWOW64 directory on the target system.

[Windows] is the Windows directory - usually C:\WINNT or C:\WINDOWS

This will improve the performance of your application on the target system.

## Deploying your 64 bit application

The current version of the library requires when deploying 64 bit applications, the Intel IPP DLLs to be deployed as well.

The DLLs are under the [install path]\LabPacks\IppDLL\Win64 directory( [install path] is the location where the library was installed).

To deploy IPP, copy the files to the [Windows]\System32 directory on the target system. [Windows] is the Windows directory - usually C:\WINNT or C:\WINDOWS  
This will improve the performance of your application on the target system.